

# 5500 W Series Genetic Analyzers Welcome Kit

## Product Insert

Publication Part Number 4477199 Revision Date 27 July 2012 (Rev. 01)

This document contains information about the iPad® and the 5500W USB drive included in the Welcome Kit. Use your iPad® to view all the 5500 W Series Genetic Analysis System documents and tutorials.

**Note:** Contact your Life Technologies Sales Representative if you are not able to access the user documents.



**WARNING! PHYSICAL INJURY HAZARD.** Do not attempt to lift or move the crated items without professional assistance. Most of the crated items are heavy. Any incorrect lifting or moving of the crated items can cause serious injury. The Life Technologies Field Service Engineers (FSEs) unpack and arrange all the components. Refer to the *5500 W Series Genetic Analyzers Site Preparation Guide* (Pub no. 4477194) for details.

---

**IMPORTANT!** Prior to operating the instrument and/or for safety and biohazard guidelines, refer to the “Safety” section in the *5500 W Genetic Analysis System User Guide* (Pub no. 4477195). For all chemicals, read the Safety Data Sheets (SDSs) and follow the handling instructions. Wear appropriate protective eyewear, clothing, and gloves.

---

## Download iTunes® on your computer

1. Go to [www.apple.com/itunes/](http://www.apple.com/itunes/).
2. Click **Download** to begin installing the latest iTunes® application.
3. Follow the prompts until iTunes® is installed.
  - Ensure that you deselect the promotional options and leave the Email Address field blank.
  - Accept the default settings.
  - Accept the license agreements, if you agree.

## Create an iTunes® account

**Note:** Before you proceed to create an iTunes® account, prepare a valid credit card (or establish a PayPal® account), a billing address, and a phone number.

1. From the menu bar click **Store**, then **Create Account**.
2. Follow the prompts until you complete creating your account.

## Synchronize your iPad® with iTunes®

1. Connect the iPad® to your computer using the supplied white iPad® USB cable.
2. If prompted, follow the prompts to register the device with Apple.
3. If prompted, to name the device, we recommend something similar to the following: **John Smith's iPad**.  
**Note:** Do not register for "MobileMe," if prompted.
4. Click **Apply** or **Sync** and allow the devices to synchronize with one another.  
**Note:** This may take a few minutes, and you may see a prompt that the device is being backed-up; this is normal.
5. If prompted to enter a password for the encrypted backup, enter the following: **1234**.
6. If not selected, select your iPad® in iTunes® under Devices.
7. Click the **Information** tab:
  - a. Scroll down and click on **Other**.
  - b. Select **Sync bookmarks with**.
  - c. Select your preferred browser.  
**Note:** If you use the Notes feature for Outlook, and if you wish to synchronize Outlook Notes via the USB cable, check the appropriate option.
8. Click the **Applications** tab and select **Sync Apps**.
9. Click the **iTunes U** tab and select **Sync iTunes U**.
10. Click **Sync**, at the bottom right-hand corner, to finalize the configuration.
11. Proceed to install **iBooks**.

---

**IMPORTANT!** In order to view or print any of the 5500 W Series Genetic Analysis System documents on your iPad®, you must first install the iBooks application on your iPad®.

---

### Install iBooks


1. Go to [www.apple.com/apps/ibooks/](http://www.apple.com/apps/ibooks/).
2. Click **Download** to begin installing the latest iBooks application.
3. Follow the prompts until iBooks is installed.

## Transfer the contents of the USB drive to iTunes®

1. Create a folder on your computer desktop and name the folder: **5500W USB**.
2. Connect the USB drive to your computer and browse to find the contents.
3. Select all of the contents of the USB drive, copy, and then paste the items into the 5500W USB folder on your computer desktop.
4. Connect your iPad® to your computer using the white iPad® USB cable.

Note: If the iTunes® application does not open automatically, open iTunes®.

5. Open the **5500W USB** folder.
6. Ensure that the 5500W USB folder and the iTunes® applications are visible on your computer screen.
7. Select the **Training Videos** folder from the 5500W USB folder.
8. Drag and drop the Training Videos folder onto the Library section of the iTunes® application.
9. Repeat steps 1 to 8 for the following folders:
  - User Documents
  - Welcome
10. Synchronize your iPad®.

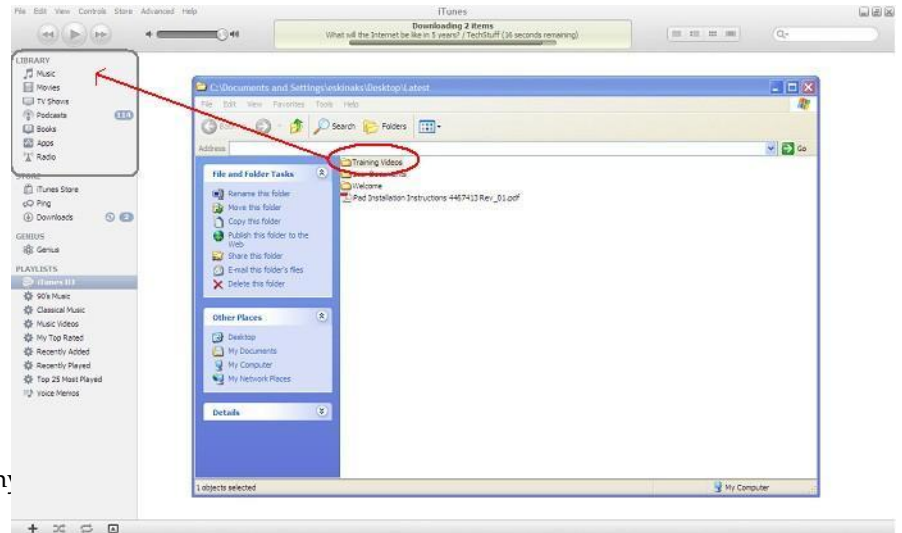

Your iPad® is now ready. Proceed to view or print any

## Access the documents

There are several methods available for accessing and using the documents. You can:

- Download the PDFs onto your computer and print the hardcopies
- Use your iPad® to review, browse, take notes, or mark (annotate) each PDF

Note: You need to synchronize the content of the USB drive with the iTunes®. See “Synchronize your iPad® with iTunes®” on page 2.

- Access the documents via the genetic analyzer [User Hub](#)

Note: To use the User Hub you must register.

## Recommended and optional additional applications for your iPad®

The following recommended and optional applications are available for *purchase*, only:

### Recommended

- GoodReader
- iThoughts
- iAnnotatePDF

### Optional

- Google Mobile App
- Adobe Ideas

For Research Use Only: Not for use in diagnostic procedures.

The information in this guide is subject to change without notice.

NOTICE TO PURCHASER: PLEASE REFER TO THE *5500 W Series Genetic Analysis Systems Limited Use Label Licenses* (Part no. 4477201) FOR LIMITED LABEL LICENSE OR DISCLAIMER INFORMATION.

#### Disclaimers

LIFE TECHNOLOGIES CORPORATION AND/OR ITS AFFILIATE(S) DISCLAIM ALL WARRANTIES WITH RESPECT TO THIS DOCUMENT, EXPRESSED OR IMPLIED, INCLUDING BUT NOT LIMITED TO THOSE OF MERCHANTABILITY, FITNESS FOR A PARTICULAR PURPOSE, OR NON-INFRINGEMENT. TO THE EXTENT ALLOWED BY LAW, IN NO EVENT SHALL LIFE TECHNOLOGIES AND/OR ITS AFFILIATE(S) BE LIABLE, WHETHER IN CONTRACT, TORT, WARRANTY, OR UNDER ANY STATUTE OR ON ANY OTHER BASIS FOR SPECIAL, INCIDENTAL, INDIRECT, PUNITIVE, MULTIPLE OR CONSEQUENTIAL DAMAGES IN CONNECTION WITH OR ARISING FROM THIS DOCUMENT, INCLUDING BUT NOT LIMITED TO THE USE THEREOF.

© 2012 Life Technologies Corporation. All rights reserved. The trademarks mentioned herein are the property of Life Technologies Corporation and/or its affiliate(s) or their respective owners. iPad and iTunes are registered trademarks of Apple, Inc. PayPal is a registered trademark of PayPal, Inc.